

CONTACT US



AJJ

HEALTHCARE MANAGEMENT PTE LTD

The related products are for display reference only, and the specific products shall prevail. Product parameters, functions, appearance and other parts of the continuous optimization, may be changed without prior notice, without further notice.

MAKE BATHING EASIER

CareX ErgoClean

Portable Patient Hygiene System



OUR PRODUCT



CareX Ergoclean

An innovative intelligent product that truly realizes portable bathing. Supports dual modes of shower and scrub, with innovative HydroSuck™ scrub, skin cleansing, and synchronized suction. Elders do not need to enter water, eliminating the hassle of movement, and can be operated by a single person, making it easy to achieve horizontal assisted bathing.

Making Bathing Simple
Make assisted bathing easier



Home care

Lack of specialization in home care



Nurse care

steps are complex, time-consuming and laborious



Elderly wish

Elders hope to take a comfortable and dignified bath



Industry status

High time and labor costs, low efficiency



1+1

ONE MAN & ONE MACHINE

CareX ErgoClean

- ✓ SCRUB
- ✓ SHOWER
- ✓ MASSAGE
- ✓ BACK-UP
- ✓ THERMOSTAT
- ✓ PORTABLE

New concept of assisted bathing

HydroSuck™: Innovative Suck Back bath
 New Bathing Concept / Firmly adhere to the skin / Synchronous suction
 Water flow control / Deep cleaning / Visible cleanliness



Instant Heating Constant Temperature

Instant heating / No need to prepare water in advance
 Save time and effort



Rotary Massage Cleaning

Skincare sponge / Strong water absorption
 Soft touch



Integrated Design

4.6L Dual bucket integrated body design
 Clear water bucket / Sewage bucket / Easy to clean



Guarding Fragile Areas of The Body

Wound avoidance / Flexible care
 Progressive Bathing / Avoid wetting the affected area


PRODUCT DETAILS



Product information

NAME: CareX ErgoClean	BUCKET VOLUME: ABOUT 4.6L
MODEL: DEN-360H	WATER TEMPERATURE REGULATION: 35°C~50°C
RATED VOLTAGE: AC 220V±22V	SPONGE ROTATION SPEED: 50RPM
RATED POWER: 1000W	COMPLIANCE WITH STANDARDS: GB4706.1-2005

CareX ErgoClean simplified control panel

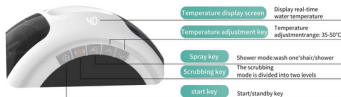

Double mode shower : Scrub / shower

Multifunctional head: supports multiple modes such as scrubbing/showering/massage.

Rotary massage function switch: rotary scrubbing for deep cleaning
Press to control water outlet and suction back

Flexible nozzle: used in shower mode

Skin friendly sponge: used in scrub mode

Shower head


Massage key

Water outlet key

Silicone head

Sponge head

Double bucket design for clean water / sewage

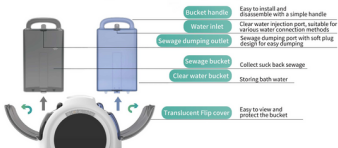

Product advantage

CareX ErgoClean Portable Patient Hygiene System

- ✓ No need to move
can be washed even in bed
- ✓ Bathe without entering the water
Suck Back bath watertight
- ✓ Partial Progressive Bathing
Protection of privacy
- ✓ Save time and effort
Bathing can be done by one person
- ✓ Instant heating
No need to prepare water in advance
- ✓ Constant water temperature
Water conservation and eco-friendly
- ✓ Convenience and Flexibility
easy use

Make Bathing Easier

VS

Traditional bath

- People who need to be moved to bed
complicated operation
- Needs to be in the water.
Excessive consumables
- Full body nudity
Be easily embarrassed susceptible to cold
- Multipayer cooperation
Increased labor consumption
- Need to prepare the venue
in advance/hot water
- High water consumption
High energy consumption
not environmentally friendly
- Scene fixation / Cumbersome
operational processes manpower
consumption / Use multiple tools
High cost of equipment maintenance

Service advantage



- 1 Improve efficiency**
one person, one machine
- 2 Comfortable Bath**
Living Water Massage Bath
- 3 Raising the bar**
Improve the standard of
assisted bathing services
- 4 worry-free and convenient**
Easier and more convenient

Door-to-door bathing aid

1、Pre-bath evaluation
Comprehensive assessment of physical condition
Skin condition consideration local cleaning

2、Preparation before bathing
Check the physical condition of the elderly
Consider environmental factors
Supplies preparation

3、Care of bath-assisting process
Consultation/ Progressive scrubbing
Standard process/ Privacy protection

4、After bathing
Auxiliary cleaning / clad /skin care
Body wetting before glazing
Check your physical condition



Make Bathing Easier

Innovation: A new way to take a convenient bath Innovative HydroSuck™ Suction bath

Efficient: Single operation is more convenient to improve the quality and efficiency of bath services

Simple: no movement, no water, more convenient, more comfortable

CERTIFICATE OF COMPETENCY

CE CERTIFICATION



SECURITY DETECTION

CE

UK
CA

RoHS

MADE IN CHINA

Independent research and development

Intelligent creation